



Vacuum Pump - Oil Free LMVP10

SALIENT FEATURES

- ✦ Electromagnetic motor and a diaphragm sealed chamber to ensure one-way gas flow
- ✦ Oil Free Pump design, relative isolation of gas circulation in the chamber
- ✦ No pollution

Applications : It can be used with corresponding equipment for aerator, nitrogen blowing instrument, solid phase extractor, thin-layer color spray, solvent filter, atmospheric sampler and other vacuum or positive pressure requirements.




E - Catalogue

HSN Code : 84141090

SPECIFICATIONS

MODEL	LMVP10
Make	LABMAN
Flow Rate	11L/min
Vacuum Pressure (Max.)	≥ 0.075Mpa
Diaphragm sealed chamber	Yes
Operating Temperature	10°C ~ 40°C
Power Wattage	55W
Power supply	AC 220V / 50Hz
Packing Dimension (LxWxH)	290 x 215 x 275mm
Gross Weight	4kg

★ Specifications are subject to change without notice.

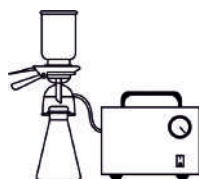

STANDARD ACCESSORIES : Main Instrument, Connecting Tube, Power Cord, Instruction Manual, Test Certificate.